

OrthoWedge **Light** Off-loading Shoe

Forefoot off-loading after surgery to correct hallux valgus or hammer toes and for healing of forefoot ulcerations



FOR RAPID HEALING AND REGENERATION

- > Light weight
- > Compression EVA outsole:
Abrasion-resistant with excellent damping properties
- > Inner sole adapted to the well-established OrthoWedge size model
- > Easy replacement of the insole systems

COMPATIBLE WITH:



PegAssist® Insole
Off-loading Insole



TwinShoe®
Balance Shoe



DARCO

The OrthoWedge Light is the optimised version of the OrthoWedge.



wide hook and loop fastener for use in a medial and lateral direction



replaceable insole



breathable mesh material



EVA outsole is easy to modify (sandable and millable)

FEATURES

- > plus for comfort: breathable upper material and removable inner sole
- > light compression EVA outsole combines shock absorbing properties and stability
- > design of the off-loading wedge (15°) produces pressure redistribution to the midfoot and rearfoot
- > the wide opening and adjustable straps ensure comfortable entry
- > wide hook and loop fastener for use in a medial and lateral direction for additional rearfoot stability
- > easy modification for targeted forefoot off-loading on the right or left foot

INDICATIONS

- > Postoperative use after surgery to correct hallux valgus or hammer/claw toes
- > Healing of forefoot ulcerations



special off-loading of the forefoot

MODES OF ACTION

- > Pressure redistribution to the midfoot or rearfoot through special off-loading and protection of the forefoot
- > Compression EVA outsole: Abrasion-resistant with excellent damping properties

COMPATIBLE WITH

SIZE ART. NO.

PegAssist® Insole

Off-loading Insole


XS	PAS0
S	PAS1
M	PAS2
L	PAS3
XL	PAS4



TwinShoe®
Balance Shoe



SIZES AND MEASUREMENTS

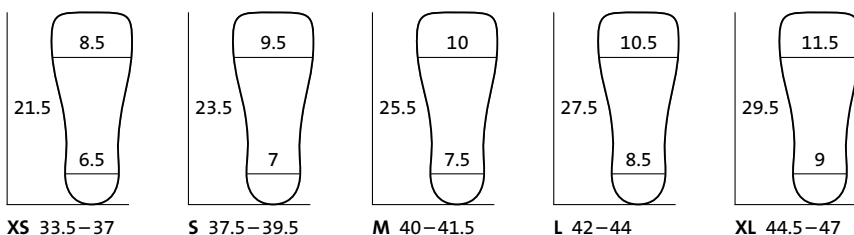

SIZE ART. NO.

OrthoWedge Light

Colour: Black

XS	OW0B-STL
S	OW1B-STL
M	OW2B-STL
L	OW3B-STL
XL	OW4B-STL